



MELAO YA TŠHOMIŠO

O swanetše o latele ditaelo tsa teko gabotse gore o hwetše dipoelo tša nnete. Ga se wa swanela go ja goba go nwa metsotse ye 15 pele o ka šomisa teko goba go šomisa didrišwa tša go hlatswa legano metsotso ye 30 pele o ka thoma teko.

TEMOŠO: Ge o le kalafong ya HIV (ARVs) o ka hwetša dipoelo tša maaka.

SEPEDI



www.oraquickhivselftest.com

ORAQUICK®
TEKO KA BOWENA YA HIV

O KA ŠOMIŠA BJANG SETLABAKELO SA ORAQUICK® HIV SELF-TEST



1

O SWANETŠE GO HWETŠA TSELA YA GO BALELA TEKO



2

Sephuthelwana se na le: **ditlabakelo tša teko, bobeele bja teko le melao ya tšhomišo.**



3

Ditlabakelo tša teko di na le sephuthelwana tše pedi.



4

Kgeila o bule sephuthelwana seo se nago le **tšhupu.**



5

Tloša sekhurumelo.




6

O se ke wa tšhola seela. O se ke wa nwa.




7

Iša tšhupu godimo ga sebeelo.

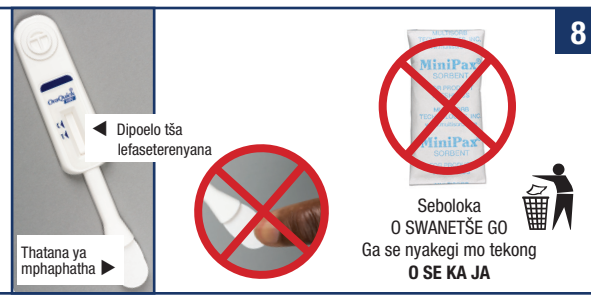


8

Kgeila o bule sephuthelwana seo se nago le sedirišwa sa teko o se tloše. **O se ke wa tshwara lefelo le la mphaphatha ka menwana ya gago. O se ke wa ja goba go nwa sebolokelo.**

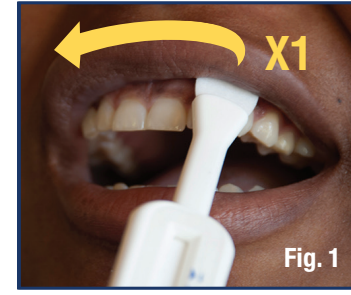

X1

Fig. 1

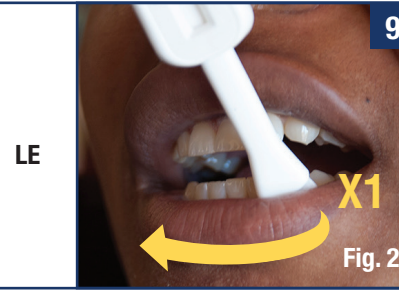

X1


Fig. 2

Kgotliša phepe ya mphaphathi go marinini a gago a ka godimo gomme o gohle gatee (fig.1) ka morago o dire se se swanago go marinini a ka fase ga tee (fig.2).



10

Tsentšha lefelo la mphatha o le sobeletše ka gare ga tšhupu go fihlela le kgoma fase.

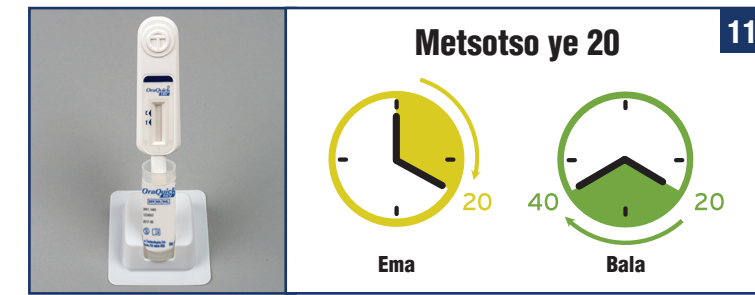

11


O e tlogelo metsotso ye 20 pele o ka bala dipoelo. O seke wa bala dipoelo ka morago ga metsotso ye 40.

GO HLATHOLLA DIPOELO



Bala dipoelo ka lefelong la go bonagala

DIPOELO TŠA GO BA PHOSETHIFI

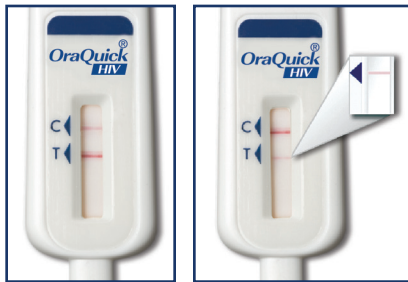

Methalo ye mebedi ya go felelela, le ge methalo ya gona ele e sa bonagale ga botse e laetša gore gona le kgonagalo ya gore o HIV phosethifi le gore o nyaka diteko tša tlaleletšo.



Ka potlako ...

Etela lefelo la go dira diteko goba tša maphelo

DIPOELO TŠA GO SE BE TŠONA



Ge go sena mothalo kgauswi le "C" (le ge go na le mothalo kgauswi le "T") goba bokhwibedu bja go dira gore o se kgone go bala dipoelo, teko ga e bereke, e swanetše go dira gape. **O swanetše go hwetša teko yengwe.**



Teko ga se ya bereka gabotse. Etela lefelo la go dira diteko tša HIV goba tsa Maphelo go dira diteko gape.

DIPOELO TŠA GO BA HIV NEKETHIFI

GE O KA BALA DITEKO PELE GA METSOTSO YE 20, DIPOLE E KA NO SEBE TŠONA



MOTHALO O TEE kgauswi le "C" go go se na le mothalo kgauswi le "T", dipoelo tša gago di nekethifi.

Nyaka diteko tša kgafetšakgafetša, ge e le gore o na le HI, dira diteko gape ka morago ga dikgwidi tše tharo.

GE O S NA LE NNETE YA DIPOELO

Ga o tsebe ka dipoelo goba go o kna bonnete ba dipoelo.

Etela lefelo la go dira diteko tša HIV goba tsa Maphelo go dira diteko gape.

GO LAHLA

Tloša thatana ya teko, bea sekhurumelo sa tšhupu ya go dira teko o lahle dilo ka moka ka gore ga ditlakala.

