



IFIKOMO FYA MIBOMFESHE

Mufwile ukukonka bwino bwino umutande uwa kubomfya pakutira mukwate ifiletumbukamo ifya mpofu. Mwilya nangu ukunwamo icili conse pali ba mintei ikumi limo nayasano elyo tamulati ukupima nangula uku bomfya umuti uuli onse uwa kuwamishya mukanwa pali ba mineti amakumi yatatu libe tamula yamba uku pima.

UBUSOKO: Ngacakuti muli pa bundapishi bwa HIV (ARV) kuti mwapimwa ifyabufi.



www.oraquickhivselftest.com

BEMBA

ORAQUICK® HIV SELF-TEST

IFYA KUBOMFYA ORAQUICK® UKUIPIMA AKASHISHI KA-HIV WE MWINE



1

MUFWILE UKUSANGA INSHILA YAKUMWENAMO INSHITA YA UKU KUPIMA



2

Umufuko mulesangwa: akakupimina, apakuteeka elyo ne nkonka ya mibomfeshe.



3

Umwaba ifya kupimina mwaba utuma tumba tubili.



4

Lepuleni ukwisula akali tumba umwaba ka chubu (tube).



5

Fumyeni inkupiko.



6

TAKULI kwitulula amenshi ayalimo. **TAKULI** ukunwa.




7

Keminikeni ka chubu (tube) apa kwimanikila.



8

Lepuleni umufuko umuli akakupiminako nokukafumya. **TAKULI** kwikatako iminwe kuli ako akapapatala. **TAKULI** kulya nangu ukumina ifyamukati ifyakukacingilila.

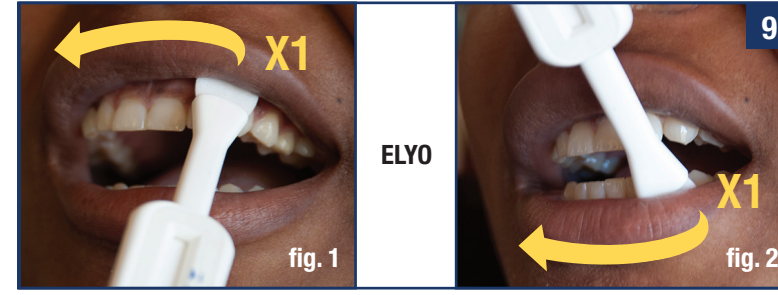

9

Shintilisheni ulubali ulwa papatala, namaka, kuya mule kululula mukati mwa ciponshi ca pa mulu, umuku umo (fig.1) nakabili mukati mwa ciponshi ca panshi, umuku umo (fig.2).



10

Bikeni ulubali ulwa papatala lonse muli ka Chubu (tube) mpaka ulubali lwapapatala lwikate panshi yaka. Chubu (tube).




11

KALEKENI UMO MWINE papite bamineti AMAKUMI YABILI elyo tamula belenga ifile tumbukamo mu mupimo. **MWI BELENGA** ifiletumbukamo mu mupimo nga capitilila pa nshita amakumi yane.

UKUBELENGA IFILE TUMBUKAMO Belengeni ifyo mwapima muncende umuli ulubuto

UKUSANGWA NAKASHISHI KA-HIV

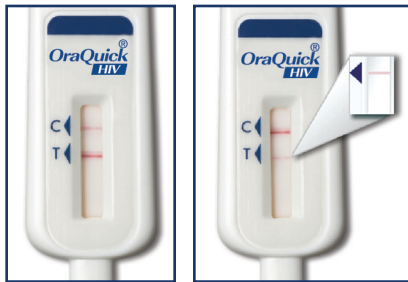

Imishilwa ibili, nangu cakuti umo taulemwenekesha, cilepilibula ukuti, kuti MWABA NA KASHISHI KA HIV elyo mulingile ukuya mukupimwa nakabili.



Mu kwa ngufyanya ...

Andalileni incende uko bapima palwa akashishi ka-HIV nangula ku kiliniki/cipatala

UKUSANGWA UWA BULA AKASHISHI KA-HIV

NGAMWA BELENGA ILYO INSHITA TAILA FIKA AMA MINITI AMAKUMI YABILI (20), NINSHI IFILETUMBUKAMO TETI FIBE IFYA CINE



UMUSHILWA UMO mupepi na "C" elyo UKWABULA UMUSHILWA mupepi na "T", capilibula ukutira tamukwete akashishi ka-HIV.

Mule pimisha libili libili. Nga cakutira mwali 'talanganapo', kapimisheni nakabili panuma/mu myeshi itatu.

IFILETUMBUKAMO IFISHILE AFWILISHA



Tapali umushilwa mupepi na "C" (nangula pali umushilwa mupepi na "T"), nangu ukutira ukukashikisha kwa nconde apo mufwile ukumwena ifya tumbukamo, umupimo taule bomba bwino elyo ulingile ukubwekeshapo. **Muyenera kutenga wina mupimo.**



Umupimo tau bombele bwino. Tandalileni/fikeni kuncende iya mupepi uko bapima palwa kashishi ka-HIV nangu ku kiliniki ukutira mupimishe nakabili.

UKUTWISHIKA IFYA TUMBUKAMO

Tamwishibe ifi tumbwikemo mu mupimo nangula muletwishika.

Tandalileni/fikeni kuncende iya mupepi uko bapima palwa kashishi ka-HIV nangu ku kiliniki ukutira mupimishe nakabili.

IFYA KUPOSA

Fumyeni aka kupiminako, bwekeshenipo inkupiko, nokuposa fyonse (mwacila bomfya), mulya muposela ifisoso lyonse.



Item# 3001-2884-70
rev. 01/18



