



LITAELO TSA TŠEBELISO EA SEHLAHISOA

U tšoaetse ho latela litaelo ka hloko hore u tle u fumane sephetho se hlakileng. U seke ua ja, ua noa letho kapa hona ho hlatsoa meno nako ea metsotso e leshome le metso e mehlano pele u etsa teko ena.

TLHOKOMELISO: Haeba u na le tšoaetso ea HIV kapa u ntse u noa li ARV u ka nna fumana sephetho se fosahetseng se bontšang kokoana-hloko e le sieo 'meleng oa hau.

SESOTHO



www.oraquickhivselftest.com

ORAQUICK[®]

HIV SELF-TEST

TSELA EA HO SEBELISA ORAQUICK[®] HIV SELF-TEST KIT



U TLA HLOKA MOKHOA OA HO SHEBA NAKO HA U NTSE U NKA TEKO/HLAHLOBO.



Sephutheloana se tšetse tse latelang: **Sesebelisoa sa thathobo, setuloana sa sesebelisoa le litaelo tsa tšebeliso ea sesebelisoa sa thathobo**



Sephutheloana sa hau se tšetse mekotlana e 'meli.



Bula mokotlana o tšetseng **botloloana**.



Bula sekoaelo sa botloloana eo.



SEKE oa qhala mokeli-keli. **SEKE** oa noa.



Kenya botlolo kahara **mokotjana** oa teng.




Bula sephutheloana se nang le **sesebelisoa sa thathobo** 'me o se ntše. **SEKE** oa tšoara sesebelisoa sa hlahlobo lehlakoreng le kenang ka hanong. **SEKE** oa ja letsoai le fumanoang ka hara sephutheloana.



Silila **bophara ba thutsoana** holima marenene a hau a kaholimo le a ka tlase hanngoe (sheba mohlala setšoantšong sa pele le sa bobeli).




Khohlelsa **bophara ba thutsoana** eo kahara botloloana ho fihlela e thetsa qetello ea botloloana.

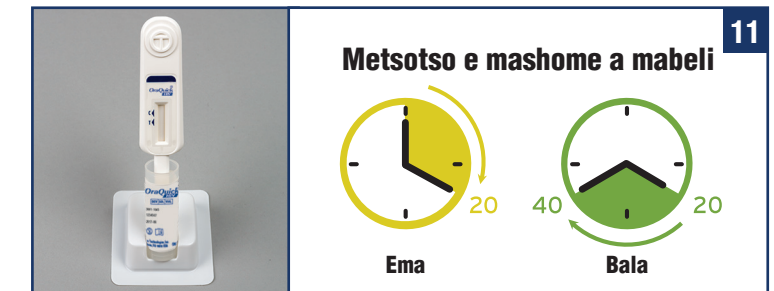

EMA METSOTSO E MASHOME A MABELI pele u ka bala sephetho. U se fetise metsotso e mashome a mane pele u ka bala sephetho.

HO TOLOKA SEPHETHO Balla sephetho khanyeng

SEPHETHO SE BONTŠANG BOTENG BA KOKOANA- HLOKO EA HIV

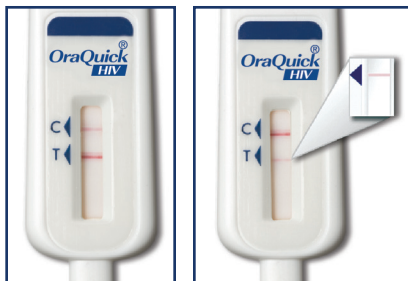

Mela e 'meli e feletseng, le haeba mola o sa teba haholo, hona ho bolela hore e kanna eaba u na le HIV 'me u lokela ho etsa hlahlobo enngoe ho netefatsa sephetho sena.



Kapele ka potlako . . .
Etela setsi se haufi sa bophelo kapa setsi sa hlahlobo tsa HIV.

SEPHETHO SE SA HLALOSEHENG



Ha mola o le sieo pela letšoao la "C" (leha o le teng pela letšoao la "T"), kapa ho hlaha bofubelu thutsoaneng eo, hona ho bolela hore sephetho ha se balehe kahoo hlahlobo e hanne 'me e tlameha ho phetoa.
U tla hloka ho etsa hlahlobo hape.



Hlahlobo ha ea etsahala hantle.
Etela setsi se haufi sa teko ea HIV kapa setsi sa bophelo 'me u etse hlahlobo hape.

SEPHETHO SE BONTŠANG BOSIEO BA KOKOANA-HLOKO EA HIV

HA U KA BALA SEPHETHO PELE METSOTSO E MASHOME A MABELI E FELA, U KA FUMANA SEPHETHO SE FOSAHETSENG



MOLA O LE MONG pela letšoao la "C" le MOLA O SIEO pela letšoao la "T", sephetho sa hau se bontša hore ha u na kokoana-hloko ea HIV.

Etsa litoko/hlahlobo khafetsa. Haeba hona le monyetla oa hore o kile oa bakotsing ea tšoaetso ea HIV, etsa litoko hape kamora khoeli tse tharo.

SEPHETHO SE SA FANENG KA BONNETE:

Ha u tsebe sephetho sa hau kapa ha se u hlakele.

Etela setsi se haufi sa hlahlobo ea HIV kapa setsi sa bophelo ho etsa hlahlobo hape.

LAHLA

Ntša thutsoana ea hlahlobo ka botloloaneng u e koahele ebe u lahla tsohle moqomong oa lithole.



E SEBELISOA HO ETSA LIPATLISISO FEELA • HA E SEBELISOE HO HLOKOMELA MOKULI

