



# ḂHO AKWANKYERE

Ḃwo se wodi nhwehwemu akwankyere no so yie pa ara sedee ebeye a wobehya nsunsuanso a eye ferenkyemm. Ansa na wobehye nhwehwemu no ase no, hwe se anye bie koraa wonnidi na wonnom biribiara sima 15 ntam ena mfa biribiara a yede hohoro anom nni dwuma sima 30 ntam ansa woahye nhwehwemu no ase.

**KOKOBO:** Se woda HIV ayarehwe (ARVs) bi so a, ebetumi aba se nsunsuanso a ebeba no nye papa.

TWI



HWE AKWANKYERE (NO)

www.oraquickhivselftest.com

# ORAQUICK<sup>®</sup>

HIV SELF-TEST

## SEDEE WODE ORAQUICK<sup>®</sup> HIV SELF-TEST KIT NO BEYE ADWUMA



**ḂWO SE WODE BERE TO NHWEHWEMU NO SO**



Boto no kura: **dee yede ye nhwehwemu no, dua a ebekura akadee no mu ne eho akwankyere.**



Dee wode beye nhwehwemu no nso kuro mmoto mmienu.



Te boto no a **tube (doroben)** no wom no.



Yi ne ti no firi so.




**NHWIE** nsuo no ngu. **NNOM.**



Fa droben no hye **dua a ebekuram** no mu.




Te boto a **nhwehwemu afidie** no wom no ano na yi firi mu. **MFA** wo nsatea nso ofa tratra no mu. **NNI** dee yede ahye no ena mmene nso.

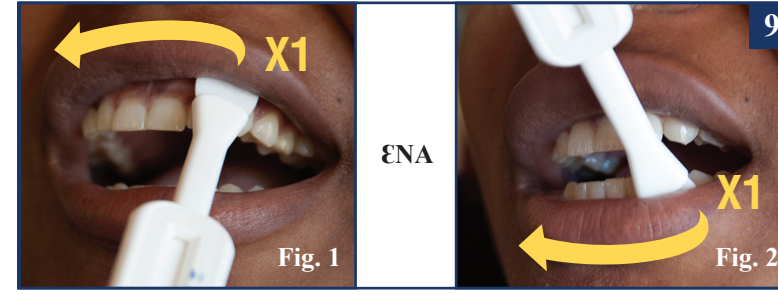

Mia Ade Tratra no dwee tare wo se akyi nam no na fa pepa wo se soro no akyi mpen baako (fig.1) ena wo se fam no akyi nam no so nso mpen baako (fig. 2).




Fa **ade tratra** no wura doroben no mu kopem se ebekopem ase.

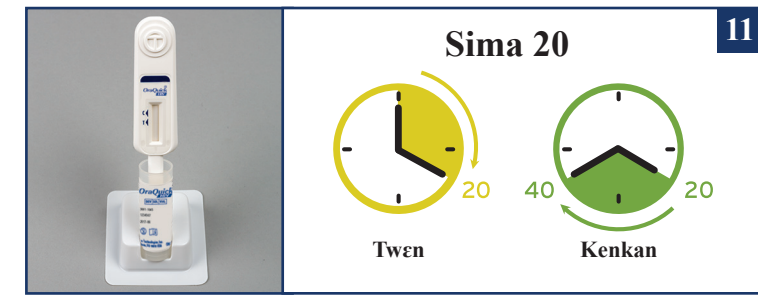

**GYAE TO HO** ma **ENI SIMA 20** ansaana wo akenkan nsunsuanso no. **NNKENKAN** nsunsuanso no se sima 40 betwam a.

## NSUNSUANSO MU NKYEREKYEREMU Kenkan nsunsuanso no wo baabi a eho ye hann pa ara

### HIV NSUNSUANSO A EKYERE SE EWOM

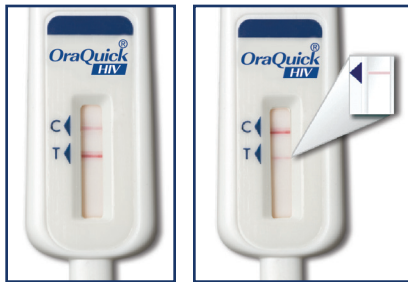

Nsensatione a ekosi mienu, se mpo se ani ye hoyaa a, ekyere se ebetumi aba se wo ANYA HIV ena ebehia se wo beko ako ye nhwehwemu bi aka ho.



*Bere a ese mu...*

**Ko HIV** Nhwehwemubea anaa Ayaresabea biara a eben wo

### NSUNSUANSO ENYE PAPA



Nsensatione biara nni "C" no ho (se mpo nsensatione bi wo "T" no ho a) anaase akyire no aye koko a ema ne kan no ye den a, nhwehwemu no nye adwuma a ewo sew aye no foforo. **Ḃwo se woma woye nhwehwemu foforo.**



Nhwehwemu no anye adwuma yie. **Ko HIV** Nhwehwemubea anaa Ayaresabea biara a eben wo ma wonye nhwehwemu no bio.

### HIV NSUNSUANSO A EKYERE SE ENNIM

## SE WOHWE ANSA SIMA 20 ABA A, EBETUMI ABA SE NSUNSUANSO NO NYE PAPA



**NSANSANEE BAAKO** wo "C" no ho a ebi nni "T" no ho a, wo HIV nsunsuanso no NNIM

Ma wantaa nye nhwehwemu. Se ebetumi aba se wowo HIV a, san so hwe wo abosome 3 ntam.

### WONTUMI NHUNU DEE NSUNSUANSO NO REPE AKYERE

Wonnim wo nsunsuanso no anaase wontumi nhunu nsunsuanso no potee.

**Ko HIV** Nhwehwemubea anaa Ayaresabea biara a eben wo ma wonye nhwehwemu no bio.

### TO GU

Yi nhwehwemu abaa no, fa ne ti no tua so nhwehwemu doroben no so na to gu sedee woto wo fi gu no ara.



**DEE YEDE YE NHWEHWEMU NKO ARA • YEMFA NHWE OYAREFOO**

